


SALTER LABS®



Salter AIRE Plus® Compressor (115V)

The Salter AIRE Plus is a durable, reliable and easy to use compressor. Aerosol delivery using the Salter AIRE Plus compressor combined with a Salter NebuTech® HDN® or Salter 8900 Series Small Volume Nebulizer provides quick, efficient and convenient delivery of aerosolized medication, minimizing medication waste and increasing patient compliance.

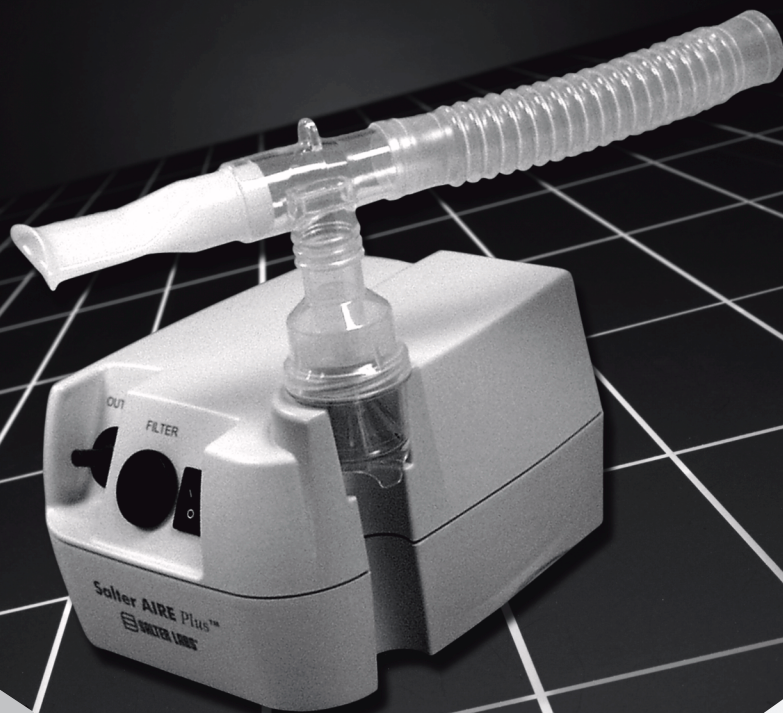

Features

- ☐ Efficient aerosol delivery
- ☐ Solid state components, 5 year warranty
- ☐ Lightweight, compact
- ☐ Piston style compressor
- ☐ Quiet, low noise level
- ☐ Matched compressor and nebulizer
- ☐ Nebulizer docking port

Benefits

- ☐ Shorter treatment time. Less medication waste
- ☐ Durable and reliable
- ☐ Easy to carry, 4.3 lbs.
- ☐ Assures strong steady pressures and flows
- ☐ Provides unobtrusive treatment
- ☐ 80% of particles in respirable range
- ☐ Added safety and convenience

Variety of uses. . . Hospital, Skilled Nursing Facility, Homecare, Clinic



slc 28


Salter AIRE Plus® Compressor (115V)

A lightweight, portable, quiet compressor for powering aerosol generating medication nebulizers. It is engineered for long term durability and cost effectiveness. The flow outlet fitting, fresh air inlet filter, and rocker style on/off switch are recessed to protect from accidental breakage. There is a convenient built in hand-grip for easy transport. A unique, tip proof nebulizer docking port is built into the case. An optional handy travel bag with extra storage for nebulizer kits is available.

Specifications

Salter AIRE Plus, Model 8350 Compressor

Size:	6.7" L x 5.7" W x 4.4" H (17 cm X 14.5 cm x 11.2 cm)
Weight:	4.3 lbs. (1.9 kg)
Electrical requirements:	115 V, 60 Hz, 1.9 A
Compressor:	Oil free cast aluminum alloy integrated uni-piston
Free Flow	10 LPM, minimum
Liter Flow	6 LPM, minimum through a Salter Labs Nebulizer
Pressure	25 PSI, minimum (312.5 KPa)
Replaceable, external filter:	White, open cell foam
UL Classified to	UL 2601-1/CAN/CSA C22.2 no.601.1
(See complete marking on product)	

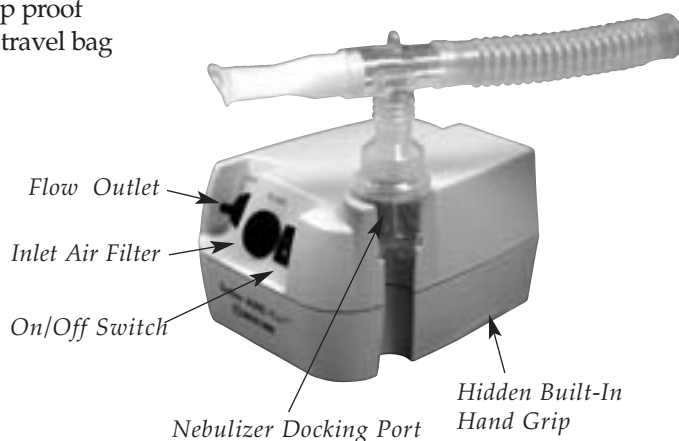

Pictured with 8900 Series Nebulizer

Ordering Information

	Part #
Compressor, Nebulizer with Tee Adapter, Mouthpiece and Reservoir Tube and 7' Tubing	8350-8900*
Same as above, plus Compressor Carrying Bag	8350B-8900*
Compressor, Nebulizer, Adult Elastic-Strap Style Aerosol Mask and 7' Tubing	8350-8924*
Compressor, Nebulizer, Pediatric Elastic-Strap Style Aerosol Mask and 7' Tubing	8350-8906*
Compressor, NebuTech® HDN® Nebulizer (Reusable) and 7' Tubing	8350-8660**
Compressor, NebuTech HDN Nebulizer (Disposable) and 7' Tubing	8350-8960**
Compressor, NebuTech HDN Nebulizer (Disposable) Adult Elastic-Strap Style Aerosol Mask and 7' Tubing	8350-8984*
Compressor, NebuTech HDN Nebulizer (Disposable) Pediatric Elastic-Strap Style Aerosol Mask and 7' Tubing	8350-8966*

Accessories	Part #
Salter AIRE Plus Compressor Carrying Bag	8258
Replacement Air Filters, 5 each per package	8501
Replacement Air Filter Cap, 1 each	8504



100 W. Sycamore Road - Arvin, CA 93203
Phone: 661-854-3166
1-800-421-0024
1-800-235-4203
Fax: 661-854-3850
Toll Free Fax: 1-800-628-4690
E-mail: salterlabs@us.salterlabs.com
www.salterlabs.com
ISO 9001 & ISO 13485

